

iSeq 100 Loading

A step-by-step guide


[iSeq System Sequencing Guide](#)




Before Loading...

- Thaw Reagent Cartridge
 - Overnight in 4°C or
 - 6 hours in room temperature water bath
- Take out flow cell from 4°C
 - At least 30 minutes before loading
- Dilute library to loading concentration (75pM or optimized concentrations)
 - Spike in 1% PhiX (default)
- Invert the cartridge 5 times then tap gently on table
- Insert flow cell into slot (B) in cartridge
- Carefully pierce foil with a pipette tip. Do not touch bottom of loading well
- Load 20µL of library into cartridge (C)

Cartridge

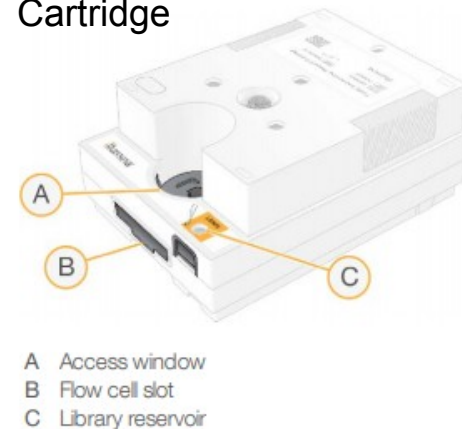

Flow cell

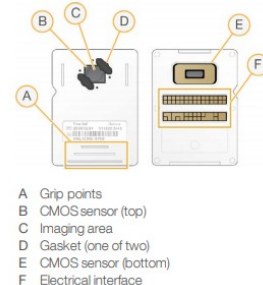



Opening Door ...



Load Cartridge?

Insert the cartridge, and then select Close Door.
Play the video for preparation instructions.



Cancel Run

Close Door



✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632



Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322

> Alerts (0)

Cancel Run

BaseSpace  SEQUENCE
HUB

michelle.tan@czbiohub.org

.....

Sign In



Configuration



Run Monitoring Only



Run Monitoring and Storage

[Skip BaseSpace Sequence Hub Sign In?](#)

✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632



Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322

> Alerts (0)

Cancel Run

Run Setup?

Run Name

Enter run name

Read Type



Single Read



Paired End

Read Cycle

151

Read 1

8

Index 1

8

Index 2

151

Read 2

> Advanced

Start Run

✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632



Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322

> Alerts (0)

Cancel Run

Run Setup ?

Run Name 20200122_Amy_K

Read Type ☐ Single Read ☒ Paired End

Read Cycle	146	12	12	146
	Read 1	Index 1	Index 2	Read 2

> Advanced

Start Run

✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632



Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322

> Alerts (0)

Cancel Run

Run Setup ?

Run Name 20200122_Amy_K

Read Type ☐ Single Read ☒ Paired End

Read Cycle	146	12	12	146
	Read 1	Index 1	Index 2	Read 2

✓ Advanced

Output Folder G:\iSeq100-01 [Browse](#)

Sample Sheet C:\Illumina\SampleSheet [Browse](#)



Start Run

≡ FS10000331 ✓(0)

✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632

Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322




> Alerts (0)

Cancel Run

Run Setup?

Run Name 20200122_Amy_K

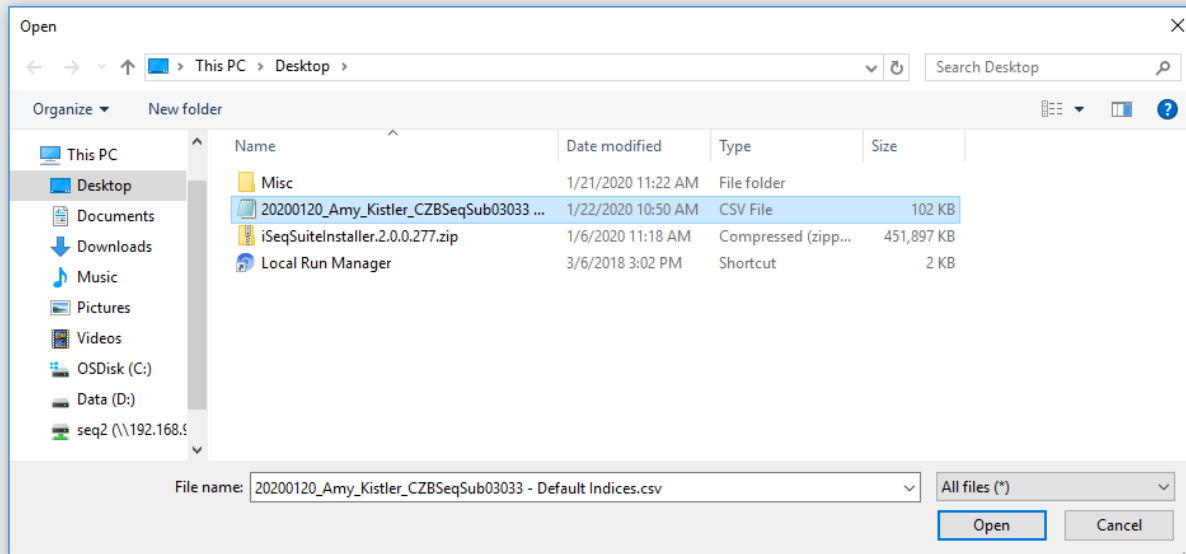

Refer to training presentation
for samplesheet instructions.

Start Run

✓ Consumables
Information

Flow Cell SN
BPA73104-1724

Lot Number
20410632



Cartridge SN
FY0026908-FC6

Lot Number
20402737

Run Cycles
322

> Alerts (0)

Cancel Run

Run Setup ?

Run Name 20200122_Amy_K

Read Type ☐ Single Read ☒ Paired End

Read Cycle	146	12	12	146
	Read 1	Index 1	Index 2	Read 2

✓ Advanced

Output Folder G:\iSeq100-01 [Browse](#)

Sample Sheet C:\Users\sbsuser.CZBIO\ [Browse](#)

Start Run




Pre-Run Checks?


The pre-run check takes about 15 minutes.
When complete, sequencing starts automatically.

20200122_Amy_K

Instrument check

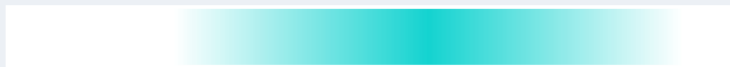

Flow check



Pre-run checks will take ~15 minutes. Run starts automatically after the checks are completed.

Cancel Run

iSeq Settings

☰ FS10000331 ✓ (0)

iSeq settings and options



← FS10000331  (0)

About

System Check

System Settings

Software Update

Process Management

Local Run Manager

Exit Application

Shut Down System



Settings iSeq™ Control Software Version 1.4.1.1700



Instrument Serial Number	FS10000331
Computer Name	ISEQ100-01
Computer IP Address	192.168.99.129
iSeq Server Version	1.4.1.654
Recipe Fragment Version	6.1.4
Total Run Count	90

Contact Us

www.illumina.com

United States
+1.800.809.ILMN
 +1.800.809.4566

Customer Service
customerservice@illumina.com

International
+1.858.202.ILMN
 +1.858.202.4566

Illumina Technical Support
techsupport@illumina.com

For license terms applicable to software included on the instrument, please see the applicable terms and conditions of sale, or terms and conditions between you and Illumina otherwise governing use of the software included on the instrument.

© 2018-2019, Illumina Inc. All Rights Reserved.

For Research Use Only

Cancel

Save

System Check

Do not remove the cartridge until the system check is complete. During mechanical test the software ejects and retracts the cartridge on



✓ Mechanical Test

✓ Thermal Test

✓ Optical Test

✓ Sensors Test

Unload

Load

Start

System Check:
Remember to put in
system test cartridge +
flow cell.

Settings

Network Access

Customization

Express Setup Settings

Illumina recommends these settings for **optimum flexibility**.



Illumina Proactive Support ?

Receive expedited support with remote diagnosis to maximize uptime.



Local Run Manager ?

Create runs, monitor status, and analyze data in a central, local location.



Remote Run Monitoring* ?

Conveniently monitor run quality and status in real time from anywhere.



Run Analysis, Collaboration, and Storage* ?

Analyze, automate, collaborate, and store data in a secure cloud space.

*Requires a free BaseSpace Sequence Hub account and internet connection.

Cancel

Set Up Manually

Use Express Settings

Using Express Settings will use **Local Run Manager** to create runs. (Different from training)

To set up runs like we do in training, please select **Set up Manually** and in the Local Run Manager setting choose **Use Manual Mode. (Next Slide)**

Settings

Network Access

Customization

Network Access

Configure internet options and other network settings.

Output Folder *

G:\iSeq100-01

[Browse](#)

Open [Network Configuration](#) to configure WiFi or Ethernet connection options

Universal Copy Service ?

Choose an account for accessing Universal Copy Service. The Windows Account setting is necessary for output folders in a directory that requires logon credentials.

☐

Local System Account

☒

Windows Account

User Name

czbiohub\sbsuser

Password

Enter password

Cancel

Save

Set output folder to write to remote server or locally. Data is saved in two places, one copy in the output folder (defined here) and another in the local drive.

File Home Share View

Pin to Quick access Copy Paste Cut Copy path Paste shortcut Move to Copy to Delete Rename New folder New item Easy access Properties Open Select all Select none Invert selection

← → > < This PC > Data (D:) > Illumina > iSeq Runs Search iSeq Runs

Quick access This PC Desktop Documents Downloads Music Pictures Videos OSDisk (C:) Data (D:) seq2 (\\192.168.99.115) (G:) Network

Name	Date modified	Type	Size
20191126_FS10000331_90_BPA73104-2820	11/27/2019 10:05 ...	File folder	
20191125_FS10000331_89_BPA73104-1718	11/26/2019 6:00 AM	File folder	
20191122_FS10000331_88_BPA73104-1119	11/23/2019 9:53 AM	File folder	
20191121_FS10000331_87_BPC29618-2020	11/22/2019 8:43 AM	File folder	
20191115_FS10000331_86_BPC29618-2508	11/16/2019 9:03 AM	File folder	
20191114_FS10000331_85_BPC29618-1310	11/15/2019 6:02 AM	File folder	
20191112_FS10000331_84_BPC29618-2211	11/13/2019 4:12 AM	File folder	
20191028_FS10000331_83_BPC29617-1734	10/29/2019 4:28 AM	File folder	
20191025_FS10000331_82_BPC29617-2528	10/26/2019 6:15 AM	File folder	
20191017_FS10000331_81_BPC29617-1030	10/18/2019 7:49 AM	File folder	
20191010_FS10000331_80_BPC29617-0727	10/11/2019 6:00 AM	File folder	
20191007_FS10000331_79_BPC29618-1821	10/8/2019 9:45 AM	File folder	
20191003_FS10000331_78_BPC29618-1109	10/4/2019 5:38 AM	File folder	
20191001_FS10000331_77_BPC29618-0820	10/2/2019 6:16 AM	File folder	
20190927_FS10000331_76_BPC29618-3119	9/28/2019 9:46 AM	File folder	
20190926_FS10000331_75_BPG80715-2135	9/27/2019 10:31 AM	File folder	
20190924_FS10000331_74_BPG80715-1737	9/25/2019 9:22 AM	File folder	
20190920_FS10000331_73_BPG80712-2237	9/20/2019 11:14 PM	File folder	
20190916_FS10000331_72_BPG80715-2537	9/17/2019 6:12 AM	File folder	
20190912_FS10000331_71_BPG80715-1331	9/13/2019 3:06 AM	File folder	
20190828_FS10000331_70_BPG80715-1835	8/29/2019 5:07 AM	File folder	
20190822_FS10000331_69_BPG80712-1637	8/23/2019 6:12 AM	File folder	
20190820_FS10000331_68_BPG80712-3227	8/21/2019 9:49 AM	File folder	
20190819_FS10000331_67_BPG80712-1334	8/20/2019 3:54 AM	File folder	
20190802_FS10000331_66_BPC29615-2336	8/3/2019 3:54 AM	File folder	

25 items | 1 item selected

11:12 AM 1/21/2020

Local copy on computer. Our sequencer also writes to our remote server.

Settings

Network Access

Customization

Instrument Customization

Name your instrument, select where to store run files locally, and other custom settings.

Instrument
Nickname

FS10000331



Audio



Off



On



Save all thumbnail images



Software Updates



Autocheck for software updates



Requires an internet connection.

Cancel

Save

← FS10000331  (0)

About

System Check

System Settings

Software Update

Process Management

Local Run Manager

Exit Application

Shut Down System



Process Management

D:\Illumina\iSeq Runs Used: 48.1 GB Available: 57.9 GB



Run Name	ID	Date	Status
20190820_Paula_S Universal Copy Service: Complete ✓	20190820_FS10000331_68_BPG80712-3227	August 20 22:50	Complete ✓ Delete Disk Space Used: 1.76 GB
20190822_Saharai_C Universal Copy Service: Complete ✓	20190822_FS10000331_69_BPG80712-1637	August 22 19:13	Complete ✓ Delete Disk Space Used: 1.78 GB
20190828_Lillian_C Universal Copy Service: Complete ✓	20190828_FS10000331_70_BPG80715-1835	August 28 18:08	Complete ✓ Delete Disk Space Used: 1.76 GB
20190912_Shay_H Universal Copy Service: Complete ✓	20190912_FS10000331_71_BPG80715-1331	September 12 19:26	Complete ✓ Delete Disk Space Used: 1.35 GB
20190916_Sreyngim_19IN109 Universal Copy Service: Complete ✓ BaseSpace Sequence Hub: Complete ✓	20190916_FS10000331_72_BPG80715-2537	September 16 19:11	Complete ✓ Delete Disk Space Used: 1.77 GB
20190920_sofonias Universal Copy Service: Complete ✓	20190920_FS10000331_73_BPG80712-2237	September 21 01:14	Complete ✓ Delete Disk Space Used: 0.25 GB
20190924_Amy_K Universal Copy Service: Complete ✓	20190924_FS10000331_74_BPG80715-1737	September 24 22:23	Complete ✓ Delete Disk Space Used: 1.69 GB

Runs deleted from Process Management are deleted from the local drive. Output folders are not deleted.

Process Management:

Can delete runs here if your computer is running out of disk space

iSeq Troubleshooting

Looking for top unknowns and re-demultiplexing

Basespace Troubleshooting → Looking for **undetermined** **reads**

>Select Analysis

OR

>Go to Projects > Select
appropriate Project > FASTQ
Generation....

SUMMARY ANALYSES BIOSAMPLES FASTQS OTHER DATASETS

Showing 1 of 1

☐	NAME
☐	FASTQ Generation 2019-12-07 01:52:16Z

basespace.illumina.com/dashboard

Apps Inbox Drive Calendar SeqSub Sub2020 Slack Quartz MyAccess AWS pacbio BaseSpace Other Bookm

Dashboard: Personal

Your FREE TRIAL has expired

Basic Tier

Usage

iCredits subscription required

Ready for more? Upgrade to a Professional subscription.

[Learn More](#)

Storage

Personal Storage in use 745.67 GB

[Details](#)

Developers

[Details](#)

Start building the next generation of BaseSpace sequencing apps.

Newsfeed

VCAT v4.0 - Support for DRAGEN and GA4GH best practices for small variant benchmarking
Author: Eric Allen, Associate Director of Bioinformatics at Illumina The...
2019-11-26

MIXCR Immune Repertoire Analyzer version 2.1.11 from MiLaboratories is now available at Illumina BaseSpace™
Youting Sun, Senior Bioinformatics Scientist at Illumina Dmitry Chudakov, CSO...
2019-11-21

BaseSpace™ CLI v1.0.0 is here!
By Swathi A. Ramani, Staff Product Manager - BaseSpace Sequence Hub If...
2019-10-07

DRAGEN™ Enrichment App - Accurate, rapid analysis for germline and somatic exome experiments
The DRAGEN Enrichment app provides...

Notifications

Share Accepted

2020-01-16 07:17
You accepted Jennifer Bohl's invitation to the Run pool 5

Share Accepted

2020-01-16 15:59
You accepted Angela Detweiler's invitation to the Run 20200116_Eric_W

Share Accepted

2020-01-14 09:35
You accepted Angela Detweiler's invitation to the Run 20200113_ReubenS

Latest Runs

[All Runs >>](#)

Running

2020-01-22 12:23
20200122_Bing_W

Running

2020-01-22 12:16
20200122_Bruce_W

Running

2020-01-22 11:10
20200122_Amy_K

Latest Analyses

[All Analyses >>](#)

Complete

2020-01-04 09:10
FASTQ Generation 2020-01-04 1 6:37:31Z

Complete

2019-12-06 18:27
FASTQ Generation 2019-12-07 0 1:52:16Z

Complete

2019-12-03 22:31
FASTQ Generation 2019-12-04 0 6:25:12Z

Unindexed Reads

Analysis: FASTQ Generation 2019-12-07 01:52:16Z

SUMMARY

REPORTS

INPUTS

FILES

General Info

Name	FASTQ Generation 2019-12-07 01:52:16Z
Application	FASTQ Generation Version: 1.0.0
Date Started	2019-12-06 17:52
Date Completed	2019-12-06 18:26
Duration	34 minutes 11 seconds
Compute Charge	0.00 iCredits
Session Type	Single Node
Size	2.58 GB
Status	Complete
Delivery	None

Logs [View Files](#)

```
Plugin assembly version information:
Common.Logging: 2.1.2.0
CsvHelper: 2.0.0.0
GenerateFASTQ: 1.1.0.64
Illumina.BaseSpace.NativeApp.Core: 0.9.0.2
Illumina.BaseSpace.NativeApp.Isis: 0.9.31.1
Illumina.BaseSpace.NativeApp.Standard: 0.9.21.1
Illumina.BaseSpace.SDK.Private: 2.0.0.45
Illumina.BaseSpace.SDK: 2.0.0.189
Illumina.TerminalVelocity: 1.0.0.141
Mono.Posix: 4.0.0.0
Newtonsoft.Json: 6.0.0.0
ServiceStack.Common: 3.9.71.0
ServiceStack.Interfaces: 1.0.0.0
ServiceStack.Text: 3.9.71.0
[APP_LOG] [2019-12-07 01:55:34.220]: Input Directory: /data/input
[APP_LOG] [2019-12-07 01:55:34.220]: Output Directory: /data/output
[APP_LOG] [2019-12-07 01:55:34.220]: Scratch Directory: /data/scratch
```

Looking at top
unknowns

Analysis: FASTQ Generation 2019-12-07 01:52:16Z

SUMMARY

REPORTS

INPUTS

FILES

NAME	CREATED	TYPE	PATH
 AdapterTrimming.txt	2019-12-06 18:24	txt	AdapterTrimming.txt
 AppSession.json	2019-12-06 18:26	json	AppSession.json
 CompletedJobInfo.xml	2019-12-06 18:24	xml	CompletedJobInfo.xml
 ConversionStats.xml	2019-12-06 18:24	xml	ConversionStats.xml
 <u>DemultiplexingStats.xml</u>	2019-12-06 18:24	xml	DemultiplexingStats.xml
 DemuxSummaryF1L1.txt	2019-12-06 18:24	txt	DemuxSummaryF1L1.txt
 FastqSummaryF1L1.txt	2019-12-06 18:24	txt	FastqSummaryF1L1.txt
 IndexMetricsOut.bin	2019-12-06 18:24	bin	IndexMetricsOut.bin
 Logging.zip	2019-12-06 18:24	zip	Logging.zip
 monitoring.log.gz	2019-12-06 18:26	gz	monitoring.log.gz
 output-20191207022414885.log	2019-12-06 18:24	log	output-20191207022414885.log
 Reports.zip	2019-12-06 18:24	zip	Reports.zip
 RunInfo.xml	2019-12-06 18:24	xml	RunInfo.xml
 SampleSheetUsed.csv	2019-12-06 18:24	csv	SampleSheetUsed.csv
 spacedock-20191207022620862.log	2019-12-06 18:26	log	spacedock-20191207022620862.log
 spacedock-infrastructure-20191207022621282.log	2019-12-06 18:26	log	spacedock-infrastructure-20191207022621282.log
 Stats.json	2019-12-06 18:24	json	Stats.json
 uploader-20191207022621949.log	2019-12-06 18:26	log	uploader-20191207022621949.log
 WorkflowError.txt	2019-12-06 18:24	txt	WorkflowError.txt
 WorkflowLog.txt	2019-12-06 18:24	txt	WorkflowLog.txt
 _metadata.json	2019-12-06 18:24	json	_metadata.json

DemuxSummaryF1L1.txt

SampleNumber	0	1	2	3	4	5	6	7	8
SampleName	None	RR043_031_stool	aR_TgLM2	RR043_032_stool	aR_CLM2	RR043_033_stool			
L1T1101	4.091762	4.524015		5.806249		4.867512		6.420163	
L1T1102	3.982026	4.492229		5.795748		4.922432		6.420995	
L1T1103	4.171473	4.538264		5.85991	4.81756	6.40814	6.240286		5.53
L1T1104	4.078548	4.532383		5.735652		4.951131		6.407804	
L1T1105	4.014317	4.520563		5.800235		4.881678		6.393103	
L1T1106	4.071557	4.484633		5.796434		4.878127		6.416164	
L1T1107	4.158956	4.532023		5.748199		4.881132		6.398271	
L1T1108	4.048977	4.499242		5.833445		4.82927	6.386219		6.26
L1T1109	4.060513	4.551185		5.874003		4.861596		6.420908	
L1T1110	4.000416	4.583814		5.850707		4.888393		6.369484	
L1T1111	3.983219	4.494364		5.825986		4.805039		6.450695	
L1T1112	4.052388	4.500664		5.789317		4.903731		6.38105	6.22
L1T1113	4.013594	4.551959		5.761779		4.878228		6.371363	
L1T1114	4.104517	4.505131		5.833227		4.892363		6.433257	
L1T1115	3.915033	4.525372		5.799429		4.89881	6.396974		6.25
L1T1116	3.98555	4.519618	5.728462		4.894692		6.391196		6.22
L1T1117	3.801323	4.593669		5.805985		4.857858		6.412584	

Close

Download (35.92 KB)

Open in text editor:

```

### Most Popular Unknown Index Sequences
### Columns: Index_Sequence Hit_Count
TCACAGAGGAGT+CACCTGTACTGA 9300
GAGTTCCTTAGCC+GGTAAGGTGAAC 3380
TCGGCGTCATCT+TCTTTCCCTACA 3120
TCACGCCAACAA+GTACAGCAGCAA 2800
CTAAACGGATCT+TCTTTCCCTACA 2540
GAGTGAGTCCTA+AACTTCGGAACA 2480
GAGTCACAAGTC+ACAGAAGGTGAA 2480
AGTGTCAAGGTC+GAACGATTGAGA 2440
AAGTCATGGCAA+CCATCTAGTGAG 2440
CAAGCTAACTGA+CGAAGATCAGCA 2400
GTGCTGTGTCA+TCCAACCAAGCA 2280
ACGACTACGTTA+TTCGGTTGGTGA 2280
CACGCCAACAAA+AGTACAGCAGCA 2160
GAGTCACAAGTC+CACAGAGGTGAA 2000
GGTTTACTATCT+TCTTTCCCTACA 1980
GTTCAGAGTCAA+TCTGCTTCACTC 1880
CAAGTCATGGCA+CATCTAGTGAGA 1880
AACCGTAAATCT+TCTTTCCCTACA 1860
GAGTGCAACCAT+TGTCACCACTAC 1820
CAAGCTAACTCG+GAAGATCAGCAA 1820
TCACGCCAACAA+AGTACAGCAGCA 1780
GAGTCACAGTCA+CACAGAAGGTGA 1780
GACTACGGTTAT+TTCGGTTGGTGA 1700
GAGTTCAGTCCA+ACGGTAAGGTGA 1680
AAGCTAACTCGA+CGAAGATCAGCA 1620
TCACAGAGGAGT+ACCTGTACTTGA 1600
GTCAGAGGTCAA+TCTGCTTCACTC 1580
AGTCACAAGTCA+CACAGAAGGTGA 1560
AGTGTCAAGTCA+CGAAGCATTGAG 1520
CAAGTCATGCAA+CCATCTAGTGAG 1500
TCACGCCAACAA+TCTTTCCCTACA 1460
AGTGTCAAGTCA+CGAAGCATTGAG 1460
TCACGCCAACAA+AGTACAGCAGCA 1420
GTTCAGAGGTCA+TCTGCTCACTCA 1400
TTCAGAGGTCAA+TCTGCTTCACTC 1380



```

Re-demultiplex with updated Sample Sheet

- You have to own the run in order to requeue analysis.
- To edit samplesheet and requeue analysis([link](#)):
 - Go to run
 - Go to samplesheet tab
 - Status (right hand side) -> Requeue -> Samplesheet
 - Edit sample sheet
 - Easier to update samplesheet in excel first, copy csv as text file and copy/paste in new samplesheet
 - Can change sample name to differentiate between two demultiplexing
 - Queue analysis
- There is a limit on how many times an analysis can be requeued. If re-demultiplexing is required again, please contact Illumina for assistance.

Examples of iSeq runs on BaseSpace

Example of a Run with Fluidic Issue ([link](#))

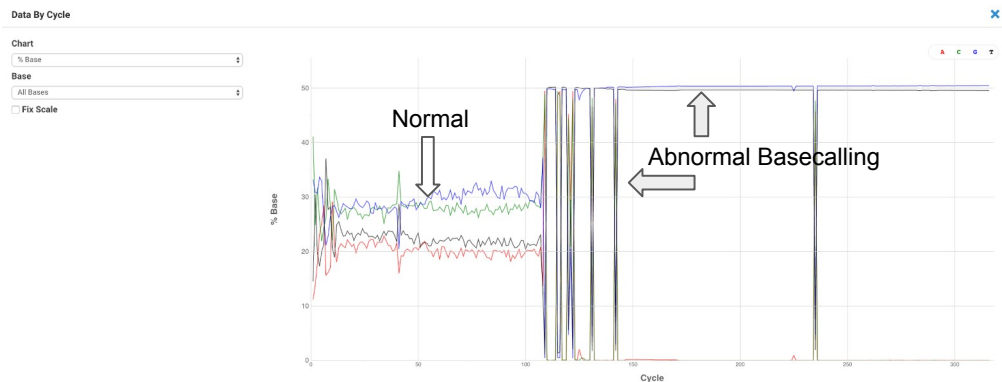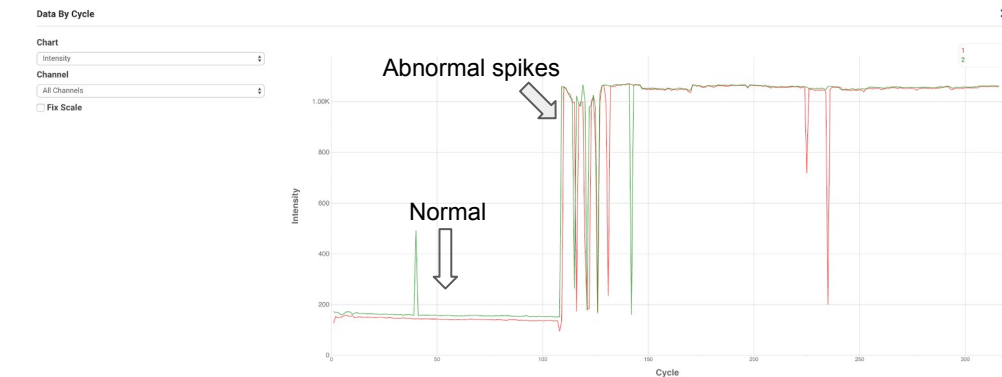

PhiX Alignment and Error Rate (usually <2%)

ALIGNED (%)	ERROR RATE (%)
2.60	21.28
0.00	0.00
0.00	0.00
0.00	0.00
1.30	21.28
1.30	21.28

In this case, the sequencer had a fluidic issue at cycle 105. The first 105 cycles have normal intensity and % base readouts, then becomes sporadic. Error rate on PhiX control is also a good indication of machine issues as error rate should be well under 2%.

Example of a Underclustered Run

% Occupied:
40-70%

Per Read Metrics

READ	CYCLES	YIELD	PROJECTED YIELD	ALIGNED (%)	ERROR RATE (%)	INTENSITY CYCLE 1	%>Q30
Read 1	135	448.27 Mbp	448.27 Mbp	0.36	1.16	88.06	90.55
Read 2 (I)	12	36.80 Mbp	36.80 Mbp	0.00	0.00	137.75	80.69
Read 3 (I)	12	36.80 Mbp	36.80 Mbp	0.00	0.00	149.19	72.08
Read 4	135	448.27 Mbp	448.27 Mbp	0.35	0.92	158.94	89.02
Non-index Reads Total	270	896.53 Mbp	896.53 Mbp	0.36	1.04	123.50	89.79
Total	294	970.13 Mbp	970.13 Mbp	0.36	1.04	133.48	88.77

20211005_MiriamS

SUMMARYSAMPLESCHARTSMETRICSINDEXING QCSAMPLE SHEETFILES

Instrument
FS10000331

88.77

Avg
%Q30

43.51

%PF

Run Status
Complete

Lane QC Status
Initial

Flow Cell Status
Initial

Latest Analysis
--

Cycles
135 | 12 | 12 | 135

Yield
970.13 Mbp

Example a Overclustered Run

% Occupied:
97.1-98.2

Per Read Metrics

READ	CYCLES	YIELD	PROJECTED YIELD	ALIGNED (%)	ERROR RATE (%)	INTENSITY CYCLE 1	%>Q30
Read 1	146	61.39 Gbp	61.39 Gbp	3.38	10.11	304.01	72.08
Read 2 (I)	12	4.66 Gbp	4.66 Gbp	0.00	0.00	1093.93	81.06
Read 3 (I)	12	4.66 Gbp	4.66 Gbp	0.00	0.00	1121.21	73.26
Read 4	146	61.38 Gbp	61.38 Gbp	1.50	4.23	207.35	71.57
Non-index Reads Total	292	122.77 Gbp	122.77 Gbp	2.44	7.98	255.68	71.82
Total	316	132.08 Gbp	132.08 Gbp	2.44	7.98	681.62	72.20

Instrument
A00111

72.20

AVG
%Q30

33.12

%PF

Run Status
Complete

Lane QC Status
Initial

Flow Cell Status
Initial

Latest Analysis
--

Cycles
146 | 12 | 12 | 146

Yield
132.08 Gbp